





Spectral Camera SWIR

OPTICAL CHARACTERISTICS	S TYPICAL SPECIFICATIONS	
Spectral range	1 000 - 2	500 nm
Spectral resolution FWHM	12 nm (30 μm slit)	
Spectral sampling	5.6 nm	
Spatial resolution	rms spot radius < 15 μm	
F/#	F/2.0	
Slith width	30 µm (50 or 80 µm optional)	
Effective slit length	9.2 mm	
ELECTRICAL CHARACTERISTI	CS	
Detector	Cryogenically cooled MCT detector	
Spatial pixels	384	
Spectral bands	288	
Pixel size	24 x 24 μm	
Detector cooling	Stirling, 25 000 h MTTF	
Optics temperature stabilization	Yes	
Camera output	16 bits CL	
SNR	1 050:1 (at max. signal level)	
Data cable	Length 5 meters	
Frame grabber	National Instruments NI 1427	
Camera control	USB / RS232	
Frame rate	450 fps (maximum full frame)	
Exposure time range	0.1 - 20 ms	
Power consumption	Nominal < 50 W	
Input voltage	Wide 24 V	
MECHANICAL CHARACTERIS	TICS	
C: (L.:\M/.:LI)	Sensor	Power supply & Control Unit
Size (L x W x H)	470 x 176 x 178 mm	300 x 190 x 130 mm
Weight	14 kg	approx. 5 kg
Body	Anodized aluminium with mounting screw holes	
Lens mount	Standard C-mount	
User adjustments	None	
Shutter	Electro-mechanical shutter for dark image acquisition	
ENVIRONMENTAL CHARACTI	ERISTICS	
Storage	-20 +50 °C	
Operating	+ 5 +40 °C, non-condensing	
SpectralDAQ support	Yes	
SDK support	Yes	
Mounting	Standard mounting see illustration, for other mounting options ask SPECIM	
Accessories	Lenses, radiometric calibration, white calibration tile, scanner stages	

ACCESSORIES

SPECIM provides various accessories for the Spectral Cameras to broaden their applicability.

Fore objective lenses, specifically designed for optimized performance in 900-2500 nm.

Lens	Focal length	FOV
OLES 15	15 mm	34 degrees
OLES 22,5	22,5 mm	23 degrees
OLES 30	30 mm	17 degrees
OLES 56	56 mm	9 degrees
OLES Macro	1:1 imaging	

Fiber optics with collection lenses or SMA connectors: from 4 to 110 input channels in one spectrometer without a moving multiplexer.

Various scanning systems: mirror scanner on rotary stage for scanning static target and outdoor scenes, and X-stage sample mover for desktop and microscope applications.


ACQUISITION SOFTWARE

SPECIM Spectral Camera SWIR is supported by Lumo software, which allows for:

- data acquisition and saving data in the hard disk
- to set camera parameters
- image visualization in real time
- to control scanner systems

Datacubes are saved in non-proprietary ENVI, Matlab and R compatible format that allows further image processing with several commercial software packages. SPECIM can also provide SDK for quick and efficient application development.

Side view

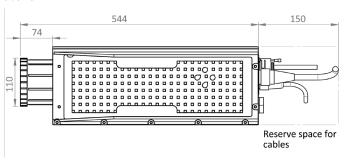

Bottom view

